


Three Steps to Start

Set up the parameter with dials, confirm with push button and start therapy with start button.



For All Patients

Adjustable AIL sensitivity.
Detect upstream occlusion.

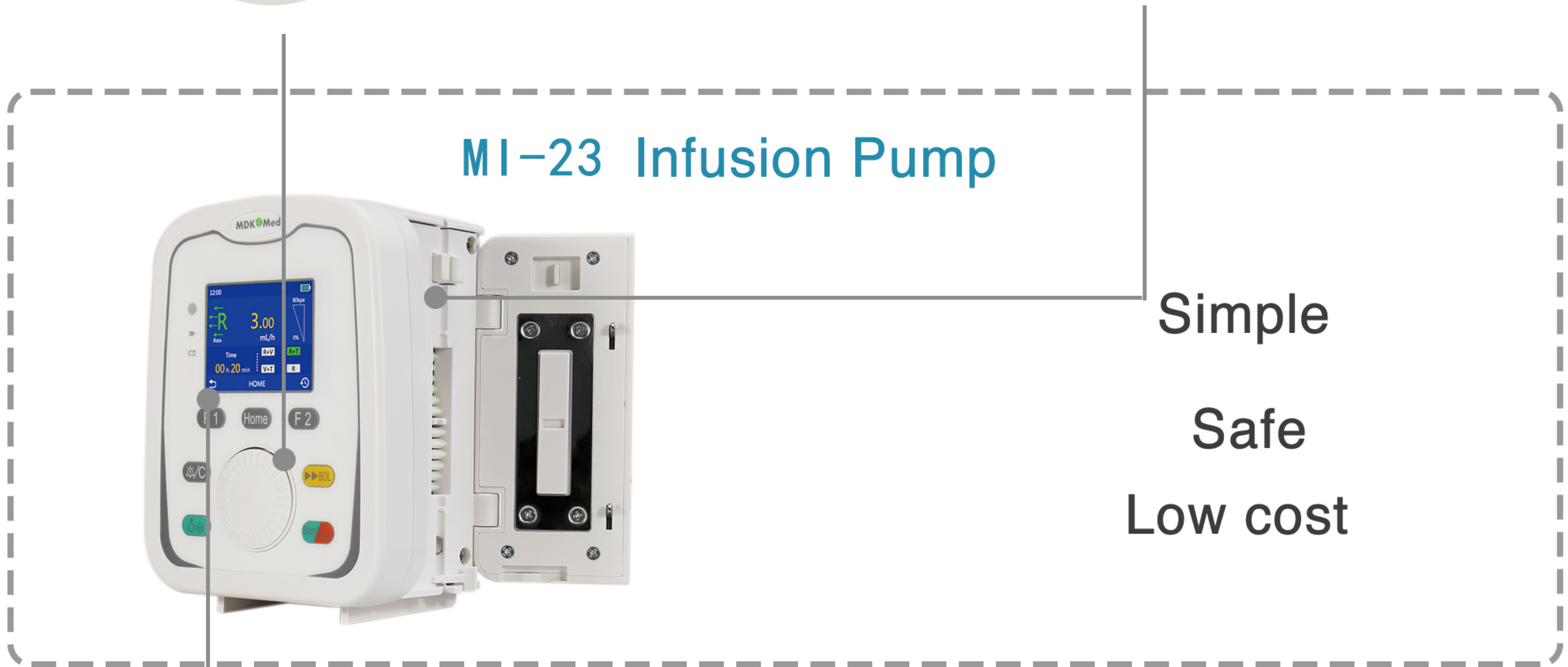

MI-23 Infusion Pump

Simple

Safe


Low cost



Touch Screen

Drug Library

Display drug name and dosage limit

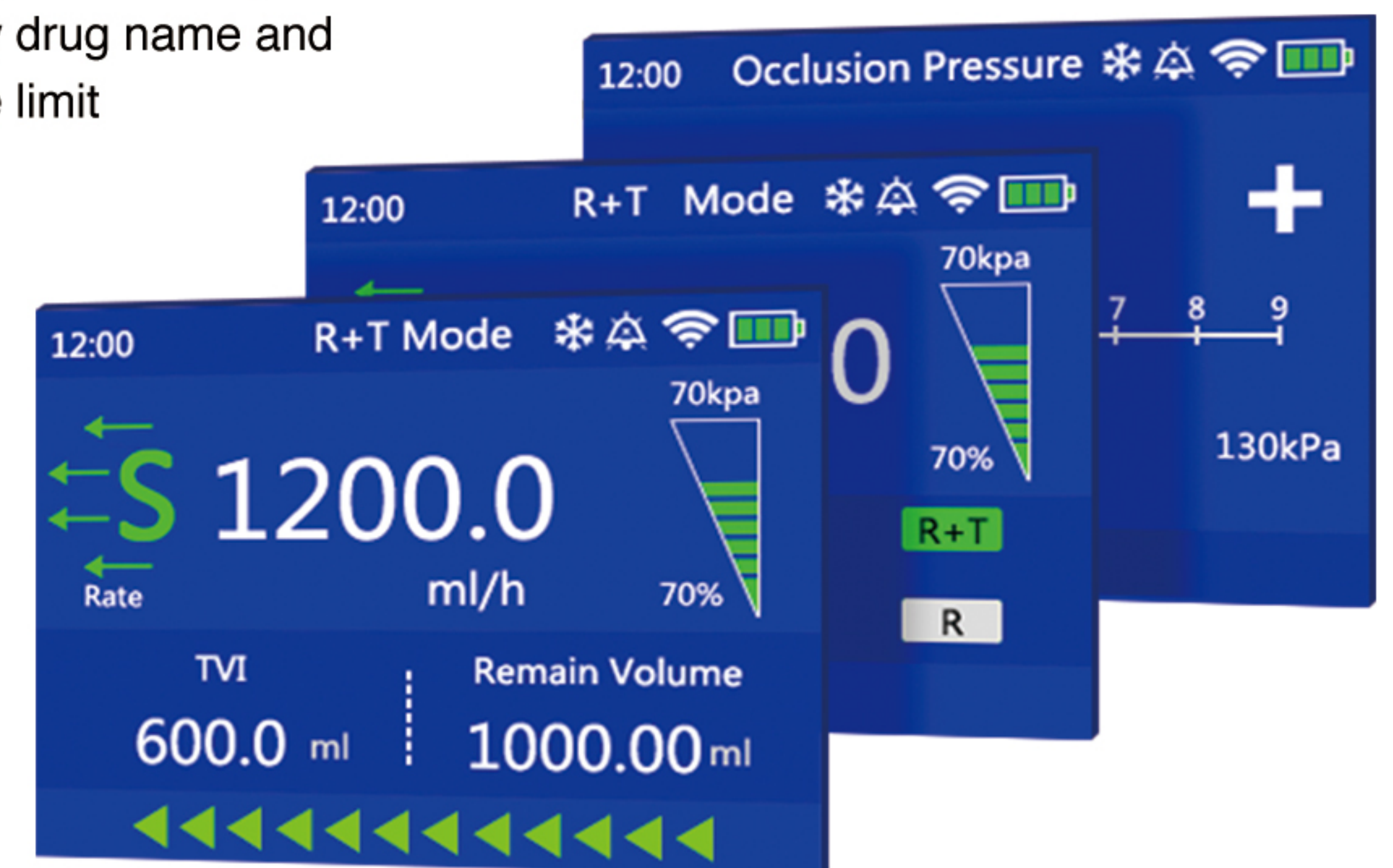

QuikMount

One click to install on the pump pole. No screw needed



Power Adapter

External power supply to reduce risk of power module overheating.



Information

Download setting and upload therapy data via WIFI



Language

Multiple languages supported



MI23 Infusion Pump

Tech Specs :

Volume Infused -----	0.01 ~ 9999.99 ml, Step by 0.01 ml
Flow Rate -----	0.01 ~ 1800 ml/h, Step by 0.01 ml/h
Accuracy -----	± 5%
Purge / Bolus Rate -----	1 ~ 1800 ml/h adjustabl, ± 20%
KVO Rate -----	Flow rate ≥ 10ml/h: KVO rate 3 ml/h Flow rate ≥1ml/h and <10ml/h: KVO rate 1 ml/h Flow rate < 1ml/h: KVO rate = Flow rate
Alarm -----	Low battery, Battery exhausted, infusion near to end alarm, Air in line, Occlusion, Door Open, Pause overtime, VTBI completion, Malfunction, A.C. and D.C. Off
Occlusion Pressure -----	Two levels: High, Low H: 100 kpa ± 30 kpa; L: 50 kpa ± 20 kpa
Power Adapter -----	A.C. input power 100 V ~ 240 V 50/60 HZ; D.C. output power 15 V/ 2A
Lithium Battery -----	Nominal Voltage: 10.8 V; Device can work over 6 hrs on battery
Infusion Pump Input Power -----	D. C. 15 V
Power -----	< 55 VA
Classification -----	Class II, Type CF, IPX2
Weight -----	1.5kg
Dimension -----	132 × 95 × 165 mm (W × D × H)

Lianying Medical Technology Co.Ltd.

Tel: 4008-808-392

Email: sales@mdkmed.com

Website: www.mdkmed.com

Company Address: Floor 20, Baoli Center, Jianggan District, Hangzhou City, Zhejiang, China.

Registered Company : Zhejiang MDKingdom Technology Co., Ltd.



Learn more